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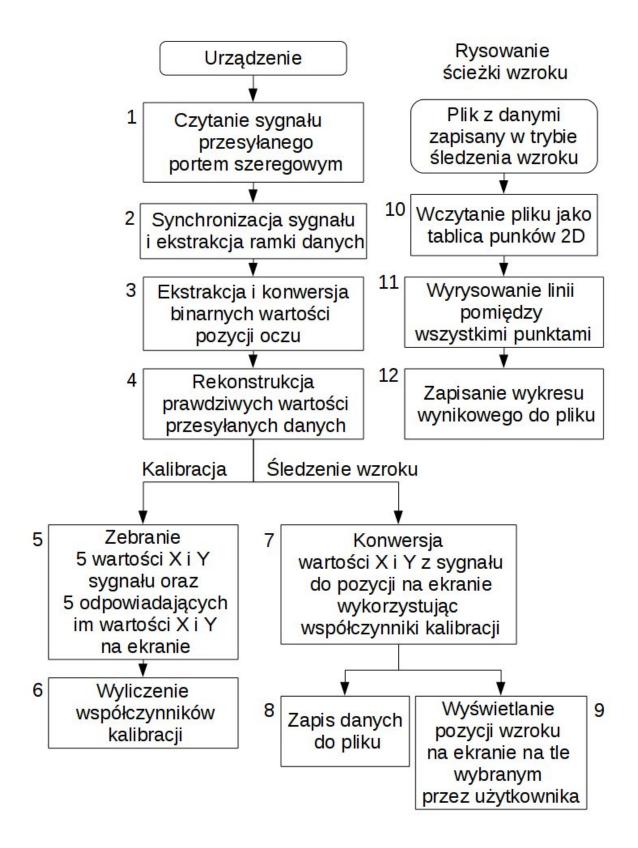
document version 15

5. Serial communication details


The JAZZ-novo system provides data to the host PC using the USB serial communication. The USB interfaces of the JAZZ-novo components utilizes the FTDI® chips. The drivers for these interfaces can be found on CD disc provided with the system.

The parameters for establishing the communication for serial transmission are following

Baud rate: 300000 bits/s

Parity: None Byte size: 8 Stop bits: 1

The JAZZ-novo system transmit data in packets using the following format.

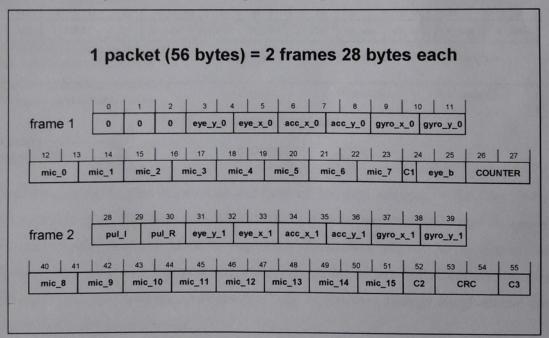Each packet contains two frames of 28 bytes each. The packets are generated with the 500 Hz frequency (2 ms per packet).


17/ June 2010

document version 15

7. Measurement window – reconstruction algorithm

Both eye position signals and pulsoxymetry signals are being measured in zoomed, movable measurement windows. Data signal samples are sent to the host PC using the 12 bits lengths format. Each sample has the offset which has to be added to the sample to obtain the sample value in the whole possible measurement range. Normally this is done automatically by the JazzRecorder software so the user does not have to take care about it. This chapter explain the reconstruction algorithm used on the PC side in case the raw data samples has to be decoded in the user's custom application.

At the start-up Jazz attempts to set the window to best-fit signal in it. After balancing, each eye position signal and pulsoxymetry signal will start in the middle of the measurement window. Movable measurement window allows tracing *Eye position* or *Pulsoxymetry* signals when they leaves measurement range. We may imagine the whole possible range of the measured singal with the sensor window moving along this range allowing to see the part of the range.

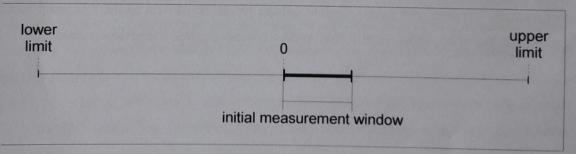
The sample values coming from the sensor to PC are fitted into the 12 bit (4096) window. We may imagine that the sensor has deducted from each sample value the offset which is actually the left window border value. Thus to obtain real sample value we have to add to it the actual value of the window left border. At the very beginning after the device is switched on we may assume that the offset is 0. Thus every sample value after adding the

17/ June 2010

offset will have the same value as before.

Let imagine that the signal value is coming close to the right border of the actual measurement window (i.e. 4095 value). When the value of the sample exceed EFF value in hexadecimal code (what is 3839 in decimal code) the next sample offset is moved to fit this sample somewhere in the measurement window. However due to the hardware implementation of the jump it is impossible to get to know the exact value of the jump in the PC software. This value has to be deducted from the following samples value. This is the moment where the reconstruction algorithm comes to play.

Basically the approach is following. Calculate hypothetical next sample value after the sample which has crossed the upper jump border EFF using the linear approximation from previous two samples (as if there was no moving window). Calculate hypothetical previous sample value to the first sample in new window using also linear approximation. Calculate first size of jump using difference between sample value from first step and sample value of the first sample in the new window. Calculate the second size of the jump using difference between last sample value in the old window and sample value from the second step. Average both jumps sizes to get the jump size closest to the real jump size and modify the offset value. As this may seem a bit complicated lets review example.

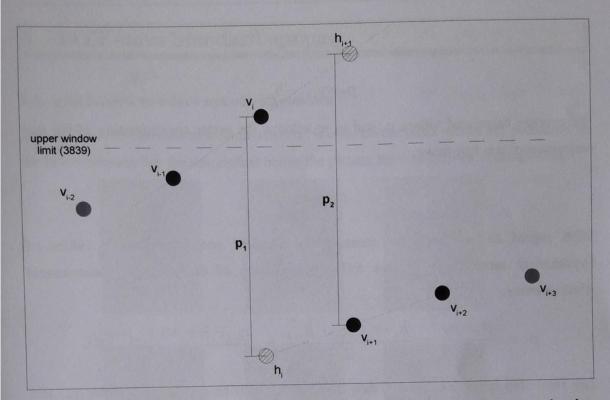As you can see on the above figure the v_i if the first value above the upper jump border. The next samples, v_{i+1} and v_{i+2} , are measured in the new, moved, measurement window (have new offset). According to our algorithm we calculate hypothetical values of sample v_{i+1} as if the measurement window had not been moved. This value is represented on the figure as sample h_{i+1} . In the similar way we calculate hypothetical value of the sample v_i in the new measurement window. This value is represented on the figure as sample h_i . Both hypothetical values are calculated using linear approximation using two consecutive samples.

$$h_{i+1} = v_i + (v_i - v_{i-1})$$

$$h_i = v_{i+1} - (v_{i+2} - v_{i+1})$$

Now we can calculate the approximation of the measurement window move based on the two pairs of samples (h_i , v_i) and (h_{i+1} , v_{i+1}).

17/June 2010

$$p_1 = h_i - v_i$$

$$p_2 = v_{i+1} - h_{i+1}$$

We average calculated values p_1 and p_2 to achieve the better approximation of the real measurement window move.

$$p = \frac{p_1 + p_2}{2}$$

With regard to equations for measurement window and equations of values of hypothetical samples we receive following equation of the averaged measurement window move.

$$p = \frac{p_1 + p_2}{2} = \frac{h_i - v_i + v_{i+1} - h_{i+1}}{2} = \frac{v_{i-1} - 3 \cdot v_i + 3 \cdot v_{i+1} - v_{i+2}}{2}$$

The final step is to calculate the new offset which will be applied to samples v_{i+1} and following. To get new offset value we have to **deduct** calculated value p from the previous offset value (let's remind here that the initial offset is 0).

$$NewOffset = OldOffset - p$$

As the move value is negative in the presented example the new offset will be greater from the old one what fits our expectation. After adding the new offset to the v_{i+1} value it will land somewhere near the h_i value. Exactly the same procedure can be used for the situation when sample value goes below the lower jump border (0FF hexadecimal or 256 decimal).

Notice! All calculations must be made using values without offset applied - just as they come from the sensor in the 12 bit (4096) range. The old offset is valid for v_i sample and all preceding, and the new offset if valid for v_{i+1} sample and all following.